


The Glia Tourniquet - Instructions For Use

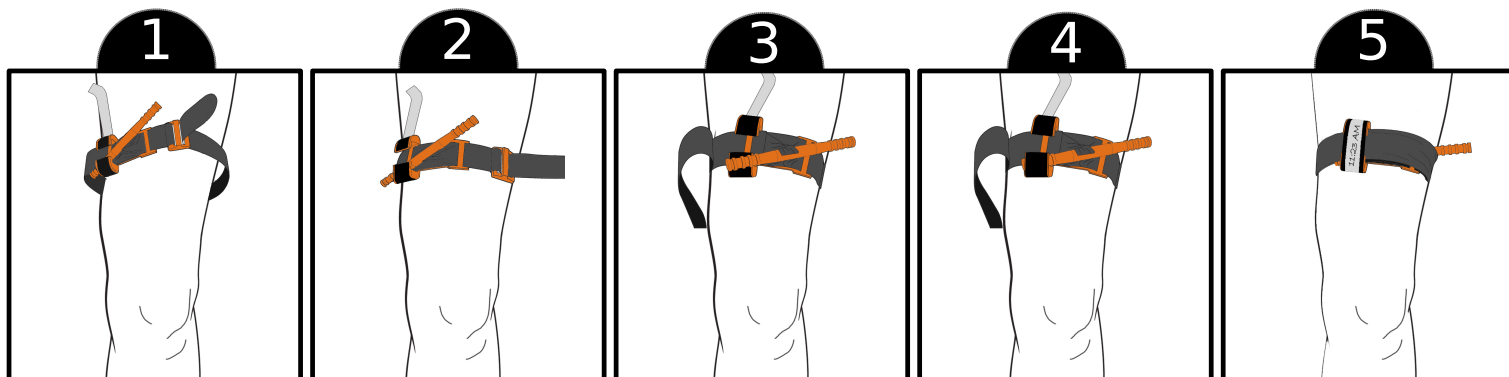

Place the injured limb through the tourniquet and place the tourniquet approximately 75mm (3") above the bleeding wound.

Where possible, the tourniquet should be placed directly on bare skin.

Tighten the strap as much as possible before securing it using the hook and loop fasteners on the strap.

Turn the windlass (rod) three half revolutions. If severe bleeding continues, secure the windlass as in the following steps and place a second tourniquet above the first as quickly as possible.

Insert the windlass inside the clip, locking it into place.

Further secure the windlass into the clip by wrapping the main strap over the windlass and securing it to itself, then wrap the white hook and loop strap over the main strap which is securing the windlass in the clip and fasten it to the black hook and loop strip on top of the clip.

Record the time of the application.


glia
equal care

Tourniquet: Model 007
54 Craig Street
London ON N6C 1E8


Tourniquet Staging

The tourniquet should come staged with the strap already threaded through the buckle.

If your tourniquet is not staged prior to use, simply thread the strap through the buckle and proceed to Step 1 above.

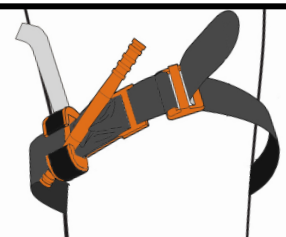

The Glia Tourniquet - Instructions For Use

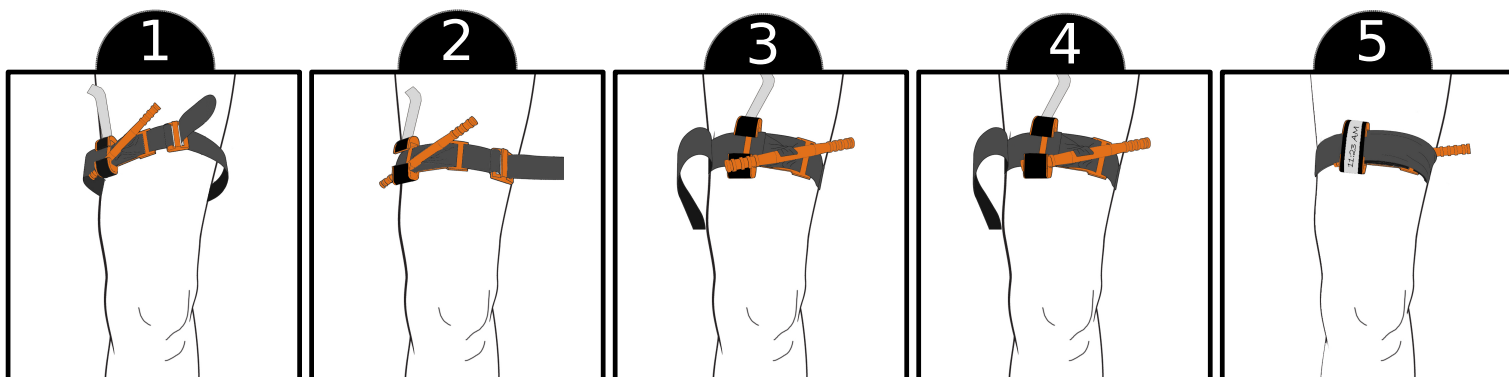

Place the injured limb through the tourniquet and place the tourniquet approximately 75mm (3") above the bleeding wound.

Where possible, the tourniquet should be placed directly on bare skin.

Tighten the strap as much as possible before securing it using the hook and loop fasteners on the strap.

Turn the windlass (rod) three half revolutions. If severe bleeding continues, secure the windlass as in the following steps and place a second tourniquet above the first as quickly as possible.

Insert the windlass inside the clip, locking it into place.

Further secure the windlass into the clip by wrapping the main strap over the windlass and securing it to itself, then wrap the white hook and loop strap over the main strap which is securing the windlass in the clip and fasten it to the black hook and loop strip on top of the clip.

Record the time of the application.

glia
equal care

Tourniquet: Model 007
54 Craig Street
London ON N6C 1E8

Tourniquet Staging

The tourniquet should come staged with the strap already threaded through the buckle.

If your tourniquet is not staged prior to use, simply thread the strap through the buckle and proceed to Step 1 above.

